



EZER™

H. pylori

Antibody Rapid Test



INTENDED USE

The EZER™ *H. pylori* Antibody Rapid Test is a lateral flow chromatographic immunoassay for the qualitative detection of antibodies to *H. pylori* in human Whole Blood/Serum/Plasma. It provides an aid in the diagnosis of infection with *H. pylori*. It's an *in vitro* diagnostic test for professional use only.



+86-571 878 18163

310018

www.genesis-ivd.com


No.139, St.10th (East), Hangzhou Economic & Technological Development Zone.

PERFORMANCE CHARACTERISTICS

Result	PCR			
	+	-	Total	
EZER™ <i>H. pylori</i> Antibody	+	142	4	146
	-	3	147	150
	Total	145	151	296


Positive Percent Agreement: 97.9%
 Negative Percent Agreement: 97.4%
 Overall Percent Agreement: 97.6%

RESULTS

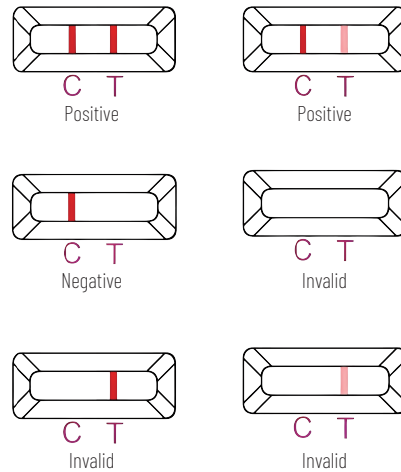


PROCEDURE

Serum or Plasma Specimen

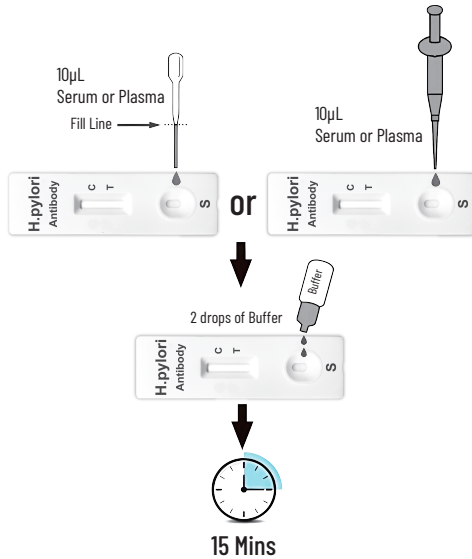

Whole Blood Specimen

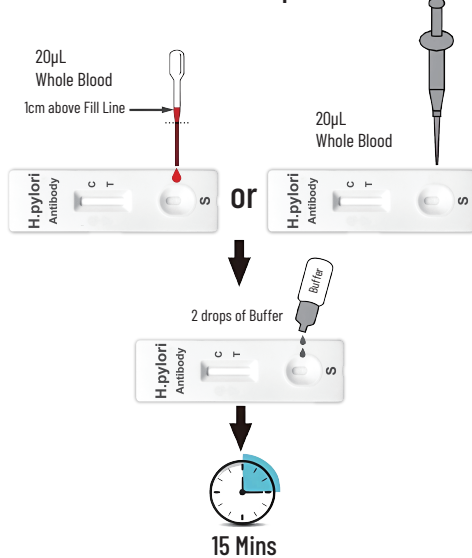

AVAILABILITY

Product	EZER™ <i>H. pylori</i> Antibody Rapid Test
Cat. No.	P221103
Test Device	40
Sample Buffer	1
Plastic Dropper	40
Package Insert	1